


# SuperCross Microcatheters

Superb Crossability in Tortuous Anatomy

# Superb Crossability in Challenging Cases

SuperCross Microcatheters provide support in tortuous anatomy and navigating bifurcated vessels.

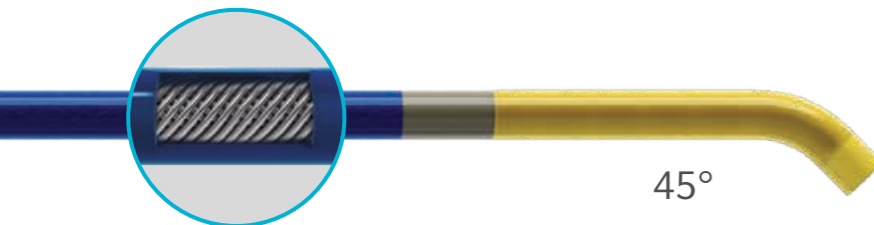


45°

## Angled Tip Catheters

### Dual Coil Shaft Construction


Provides excellent torque response, flexibility, pushability and kink resistance

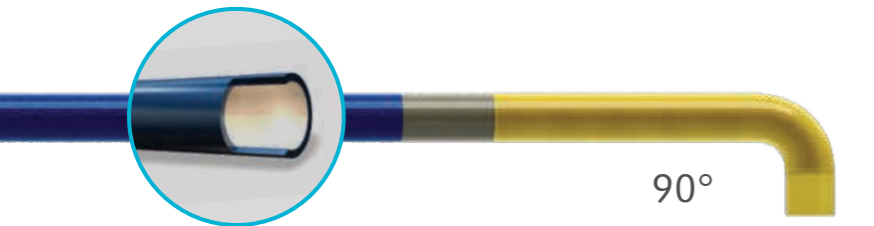

90°

### PTFE Inner Layer

Allows for outstanding guidewire movement and delivery

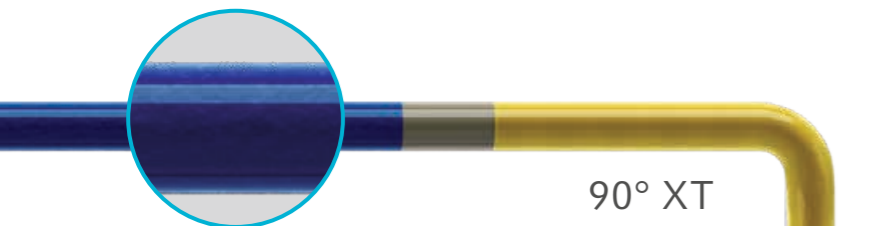

90° XT

XT=90° Extended Tip for Secure Cannulation

### Hydrophilic Coating on the Distal 80 cm

Enables exceptional trackability around tight bends and tortuous anatomy



120°

### Embedded Platinum/Tungsten Coil

Provides enhanced visibility along the entire angled tip



Supercross FT

### Stainless Steel Braid

For flexibility, pushability and kink resistance

## Flexible Tip (FT) Catheter

### Flexible Tip (FT)

For increased trackability through tortuous anatomy


## Ordering Information

### SuperCross Microcatheters

The SuperCross Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5302	0.014"	45°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5303	0.014"	45°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5304	0.014"	90°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5305	0.014"	90°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5306	0.014"	120°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5307	0.014"	120°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5308	0.014"	90° XT	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5309	0.014"	90° XT	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm

Packaged in quantities of 1 unit per box.

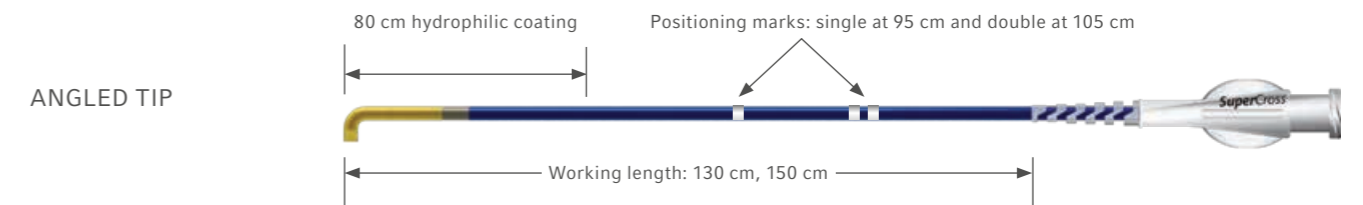


### SuperCross Microcatheter FT

The SuperCross Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5340	0.014"	FT	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5341	0.014"	FT	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm

Packaged in quantities of 1 unit per box.

FT= Flexible Tip

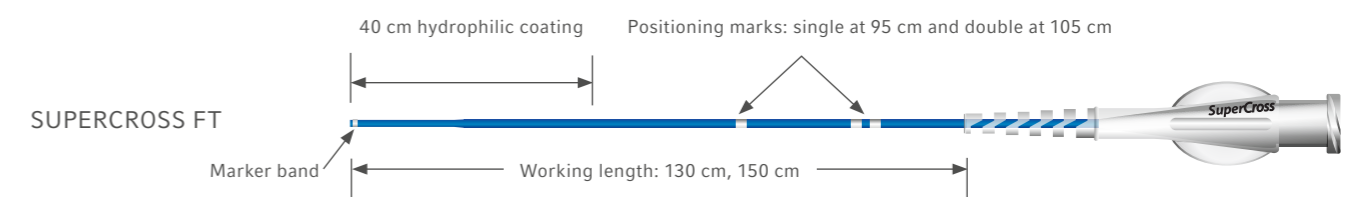

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit [teleflex.com](https://teleflex.com).

Teleflex is the home of Arrow, Deknatel, LMA, Pilling, QuikClot, Rüschi, UroLift, and Weck – trusted brands united by a common sense of purpose.

### Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

### Regional Offices

**United States:** Phone +1 919 544 8000, Toll Free 866 246 6990, [cs@teleflex.com](mailto:cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**Latin America:** Phone +1 919 433 4999, [la.cs@teleflex.com](mailto:la.cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**International:** Phone +353 (0)9 06 46 08 00, [orders.intl@teleflex.com](mailto:orders.intl@teleflex.com), Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

**Australia** 1300 360 226

**Austria** +43 (0)1 402 47 72

**Belgium** +32 (0)2 333 24 60

**Canada** +1 (0)800 387 9699

**China (Shanghai)** +86 (0)21 6060 7000

**China (Beijing)** +86 (0)10 6418 5699

**Czech Republic** +420 (0)495 759 111

**France** +33 (0)5 62 18 79 40

**Germany** +49 (0)7151 691 3004

**Greece** +30 210 67 77 717

**India** +91 (0)80 4093 4790

**Italy** +39 0362 58911

**Japan** +81 (0)3 6632 3600

**Korea** +82 2 536 7550

**Mexico** +52 55 5002 3500

**Netherlands** +31 (0)88 00 215 00

**New Zealand** 0800 601 100

**Poland** +48 22 4624032

**Portugal** +351 225 491 051

**Singapore** (SEA non-direct sales countries) +65 6439 3000

**Slovak Republic** +421 (0)3377 254 28

**South Africa** +27 (0)11 807 4887

**Spain** +34 91 198 84 31

**Switzerland** +41 (0)31 818 40 90

**United Kingdom** +44 (0)1494 53 27 61

**For more information, please visit [teleflex.com](https://teleflex.com)**

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

Teleflex, the Teleflex logo, Arrow, Deknatel, LMA, Pilling, QuikClot, Rüschi, SuperCross, UroLift and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 03/2022.

© 2023 Teleflex Incorporated. All rights reserved.

MCI-2018-0396-EN · REV 1 · 09 23 PDF