

SUPER ARROW-FLEX® SHEATH INTRODUCER

Flexibility you can see. Strength you can feel.



COMPLETE VASCULAR ACCESS WITHOUT KINKING OR COLLAPSING

THE FLEXIBLE, KINK-RESISTANT SUPER ARROW-FLEX SHEATH INTRODUCER

Super Arrow-Flex gives you complete vascular access, thanks to its superior design and construction. From its easy insertion and handling to its highly visible radiopaque tip,* Super Arrow-Flex is designed to deliver optimal performance from start to finish. With a wide range of sheath sizes and lengths, Super Arrow-Flex gives you flexibility in choices for multiple vascular applications.

MAXIMUM EFFECTIVENESS EVEN IN THE MOST CHALLENGING CASES

- Exclusive coil-wire design allows the sheath to flex at any point in any direction without kinking or losing support
- Excellent steerability helps negotiate tortuous anatomies
- Tapered dilator provides exceptional pushability to the contralateral leg and Superficial Femoral Artery (SFA)

EASY PLACEMENT AND INSERTION

- Radiopaque tip marker enhances visualization during insertion and enables precise deployment of a catheter just beyond the sheath tip*
- Highly radiopaque sheath, yet spacing between coils enables visualization through the sheath
- Smooth transition between the sheath and dilator helps to prevent peel-back
- Hydrophilic surface helps facilitate smooth insertion
- Valve lubrication eases catheter placement through the sheath


RANGE OF PRODUCTS FOR INTERVENTIONAL PROCEDURES

- Specific designs suitable for renal, carotid and transseptal procedures
- Wide range of products providing reliable access for crossover, TIPS and antegrade brachial procedures
- Variety of product options from 5 Fr. To 11 Fr. and 7.5 cm to 100 cm
- Sheath available with and without guide wire
- Color-coded for convenience

SAFETY AND PEACE OF MIND

- Ability to be used as a guiding sheath may eliminate the need for a guiding catheter (includes obturator cap)
- Color-coded hubs for easy identification of sheath sizes
- Tri-cusped hemostasis valve facilitates catheter manipulation with minimal bleed-back
- Sheath and dilator together to enable simultaneous, secure insertion

* In all 35 cm and longer except 5 Fr.

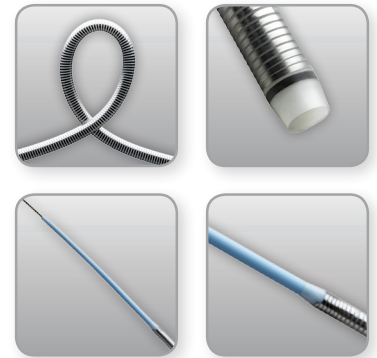

EACH CP SET INCLUDES:

- One Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One 17 3/4" (45 cm) long dual purpose "J" tip spring-wire guide in holder with tip straightener. Short-length sheath set (P) has .021 dia. and all others have .035" dia. wire guide
- One Vessel Dilator with SnapLock™ Feature
- One Obturator Cap

EACH CL SET INCLUDES:

- One Super Arrow-Flex Sheath with Integral Side Port/Hemostasis Valve and attached 3-Way Stopcock
- One Vessel Dilator with SnapLock Feature
- One Obturator Cap

**ORDERING INFORMATION****SUPER ARROW-FLEX SHEATH INTRODUCER SETS WITH WIRE GUIDES**

COLOR CODE	PRODUCT NO.	SHEATH SIZE	SHEATH LENGTH	SETS PER CASE
	CP-07511-P	5 Fr.	7.5 cm	10
	CP-07511	5 Fr.	11 cm	10
	CP-07611-P	6 Fr.	7.5 cm	10
	CP-07611	6 Fr.	11 cm	10
	CP-07711	7 Fr.	11 cm	10
	CP-07811	8 Fr.	11 cm	10
	CP-07911	9 Fr.	11 cm	10
	CP-07011	10 Fr.	11 cm	10

SUPER ARROW-FLEX SHEATH INTRODUCER SETS WITHOUT WIRE GUIDES

COLOR CODE	PRODUCT NO.	SHEATH SIZE	SHEATH LENGTH	SETS PER CASE
5 FR. INTRODUCERS				
	CL-07511	5 Fr.	11 cm	10
	CL-07524	5 Fr.	24 cm	5
	CL-07545	5 Fr.	45 cm	5
	CL-07565	5 Fr.	65 cm	5
	CL-07590	5 Fr.	90 cm	2
6 FR. INTRODUCERS				
	CL-07611	6 Fr.	11 cm	10
	CL-07624	6 Fr.	24 cm	5
	CL-07635	6 Fr.	35 cm	5
	CL-07645	6 Fr.	45 cm	5
	CR-07645*	6 Fr.	45 cm (Renal)	5
	CL-07665	6 Fr.	65 cm	5
	CL-07690	6 Fr.	90 cm	2
7 FR. INTRODUCERS				
	CL-07711	7 Fr.	11 cm	10
	CL-07724	7 Fr.	24 cm	5
	CL-07735	7 Fr.	35 cm	5
	CL-07745	7 Fr.	45 cm	5
	CR-07745*	7 Fr.	45 cm (Renal)	5
	CR-07745-NT†	7 Fr.	45 cm (Renal)	5
	CL-07765	7 Fr.	65 cm	5
	CL-07780	7 Fr.	80 cm	2
	CL-07790-R*	7 Fr.	90 cm (Carotid)	2
	CL-07700	7 Fr.	100 cm	2
8 FR. INTRODUCERS				
	CL-07811	8 Fr.	11 cm	10
	CL-07824	8 Fr.	24 cm	5
	CL-07835	8 Fr.	35 cm	5
	CL-07845	8 Fr.	45 cm	5
	CT-07860	8 Fr.	60 cm (Transseptal)	2
	CL-07865	8 Fr.	65 cm	5
	CL-07880	8 Fr.	80 cm	2
	CL-07880-R*	8 Fr.	90 cm (Carotid)	2
	CL-07800	8 Fr.	100 cm	2
9 FR. INTRODUCERS				
	CL-07911	9 Fr.	11 cm	10
	CL-07924	9 Fr.	24 cm	5
	CL-07965	9 Fr.	65 cm	5
	CL-07980	9 Fr.	80 cm	2
	CL-07900	9 Fr.	100 cm	2
10 FR. INTRODUCERS				
	CL-07011	10 Fr.	11 cm	10
	CL-07024	10 Fr.	24 cm	5
	CL-07035	10 Fr.	35 cm	5
	CL-07045	10 Fr.	45 cm	5
	CL-07065	10 Fr.	65 cm	5
	CL-07080	10 Fr.	80 cm	2
11 FR. INTRODUCERS				
	CL-71165	11 Fr.	65 cm	5
	CL-71180	11 Fr.	80 cm	2

*These products include a special Tuohy Borst "Y" Adapter with Side Arm Extension

†Tuohy-Borst "Y" Adapter is not included with this product

Caution: US Federal law limits this device to sale by or on order of a physician

Teleflex is a global provider of medical products designed to enable healthcare providers to protect against infections and improve patient and provider safety. The company specializes in products and services for vascular access, respiratory, general and regional anesthesia, cardiac care, urology and surgery. Teleflex also provides specialty products for device manufacturers.

The Teleflex family of brands includes ARROW®, BEERE MEDICAL®, DEKNATEL®, GIBECK®, HUDSON RCI®, KMEDIC®, PILLING®, PLEUR-EVAC®, RÜSCH®, SHERIDAN®, SMD®, TAUT®, TFX OEM®, VASONOVA™ and WECK®, all of which are trademarks or registered trademarks of Teleflex Incorporated or its affiliates. Teleflex Incorporated (NYSE:TFX) has annual revenues of approximately \$1.4 billion and customers in more than 130 countries.

For detailed information, see teleflex.com

TELEFLEX MEDICAL

PO Box 12600

Research Triangle Park, NC 27709

Toll Free 866.246.6990

Phone 919.544.8000

Intl 919.433.8088

Teleflex, Arrow, SnapLock and Super Arrow-Flex are trademarks or registered trademarks of Teleflex Incorporated or its affiliates.
© 2012 Teleflex Incorporated. All rights reserved. 2012-0722