



Rüsch®
CrystalClear™ Plus
 Tracheostomy Tube

Supports Effective Tracheostomy Care

The single-use Rüsch® CrystalClear™ Plus Tracheostomy Tube portfolio from Teleflex® is comprised of cuffed, uncuffed, or fenestrated (cuffed and uncuffed) tracheostomy tubes, which are sold with two single-use, flexible inner cannulas. The range is ideally suited for a new Tracheostomy stoma, ‘first tube in’ situations. Manufactured from thermo-sensitive PVC tubing, and with a thin, flexible, PVC flange, Rüsch® CrystalClear™ Plus Tracheostomy Tubes have been designed for patient comfort and effective stoma care¹. The flexibility of the flange enables it to conform to the anatomy with minimal skin contact, when a stoma dressing is not required; it also helps to minimise tube movement during care, and provides ease of access to the stoma site for cleaning, to help reduce the risk of infection².

Designed with Patient Safety in Mind

Rüsch® CrystalClear™ Plus Tracheostomy Tubes have size-specific colour-coded flexible inner cannulas that correspond to coloured writing on the flange. This helps to avoid the inadvertent use of an inappropriately sized inner cannula, which could increase ‘work of breathing’ and potentially compromise patient safety. Moreover, the flexible inner cannula can be removed and replaced, helping to prevent the accumulation of biofilm and secretions, which could occlude the outer cannula and restrict ventilation. Additional safety features include a swivel on the 15 mm connector to reduce torque on the ventilator circuit, and the 95° curvature of the tube, which facilitates a secure fit in new stomas to reduce the risk of accidental decannulation.

Benefits



The Clinician
 Colour-coded inner cannulas aid size identification, for increased peace of mind



Your Institution
 Designed to reduce the risk of inadvertent tube displacement in the intensive care unit




The Patient
 Design features assist in maintaining a patent airway

Rüsch® CrystalClear™ Plus Tracheostomy Tube

15 mm Swivel Connector
Designed to reduce torque with ventilator circuits ...

Cuff
Tubes are available both with or without a cuff ...

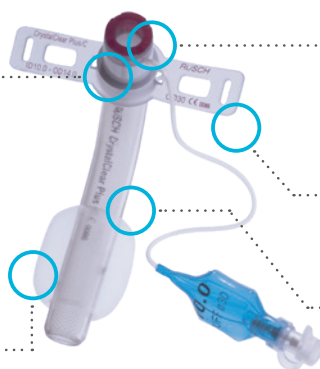

Inner Cannula
Size-specific colour coding to aid size identification

Thin, Flexible Flange
Facilitates effective stoma care

95° Shaft
Promotes a secure fit to help prevent accidental decannulation

Rüsch® CrystalClear™ Plus Tracheostomy Tube Fast Facts:

Single-use portfolio to help reduce the risk of infection associated with reusable tracheostomy tubes.

Sold with two flexible, colour-coded inner cannulas.

Flexible Flange to allow good access to stoma.

Compatible with the Rüsch® PercuTwist® Percutaneous Tracheostomy Dilator.

Rüsch® CrystalClear™ Plus Tracheostomy Tube

PRODUCT CODE*	PRODUCT CODE**	SIZE / I.D. (MM)	O.D. (MM)	LENGTH (MM)	ANGLE Ø APPROX.	CUFF Ø (MM)	INNER CANNULA	COLOUR CODE	QTY
UN-FENESTRATED									
121300-000035	121303-000035	3.5	5.3	47.0	95°	11.0	-	-	5
121300-000040	121303-000040	4.0	5.9	54.0	95°	12.0	-	-	5
121300-000045	121303-000045	4.5	6.5	54.0	95°	14.0	-	-	5
121300-000050	121303-000050	5.0	7.0	57.0	95°	14.0	-	-	5
121300-000055	121303-000055	5.5	7.5	58.0	95°	15.0	-	-	5
121300-000060	121303-000060	6.0	8.7	64.0	95°	19.0	6.0	Orange	5
121300-000065	121303-000065	6.5	8.7	64.0	95°	19.0	6.5	Maroon	5
121300-000070	121303-000070	7.0	10.0	70.0	95°	21.0	7.0	Blue	5
121300-000075	121303-000075	7.5	10.7	75.5	95°	25.0	7.5	Yellow	5
121300-000080	121303-000080	8.0	11.3	81.0	95°	28.0	8.0	Lilac	5
121300-000085	121303-000085	8.5	11.3	81.0	95°	28.0	8.5	Green	5
121300-000090	121303-000090	9.0	12.7	87.0	95°	28.0	9.0	Brown	5
121300-000095	121303-000095	9.5	12.7	87.0	95°	28.0	9.0	Navy	5
121300-000100	121303-000100	10.0	14.0	95.5	95°	30.0	10.0	Red	5

* Without Cuff. ** With Cuff. *** Available December 2018

Rüsch® CrystalClear™ PDT Set

PRODUCT CODE	SIZE / I.D. (MM)	O.D. (MM)	LENGTH (MM)	ANGLE Ø APPROX.	CUFF Ø (MM)	COLOUR CODE	QTY
UN-FENESTRATED WITH CUFF							
121304-000070	7.0	10.0	70.0	95°	21.0	Blue	1
121304-000080	8.0	11.3	81.0	95°	28.0	Lilac	1
121304-000090	9.0	12.7	87.0	95°	28.0	Brown	1

Rüsch® CrystalClear™ Accessories

PRODUCT CODE	DESCRIPTION	QTY
121311	Shower cap for Rüsch® CrystalClear™ Plus	4
121312	Cough cap for Rüsch® CrystalClear™ Plus	4
121313	Sealing cap for Rüsch® CrystalClear™ Plus	4

Rüsch® CrystalClear™ Inner Cannula

PRODUCT CODE	SIZE	I.D. (MM)	O.D. (MM)	COLOUR CODE	QTY
121310-000060	6.0	5.0	5.7	Orange	10
121310-000065	6.5	5.0	5.7	Maroon	10
121310-000070	7.0	6.0	6.7	Blue	10
121310-000075	7.5	6.5	7.3	Yellow	10
121310-000080	8.0	7.0	7.7	Lilac	10
121310-000085	8.5	7.0	7.7	Green	10
121310-000090	9.0	8.0	8.7	Brown	10
121310-000095	9.5	8.0	8.7	Navy	10
121310-000010	10.0	8.8	9.7	Red	10

References:

1. Wilkinson KA, Martin IC, Freeth H, Kelly K, Mason M. *NCEPOD: On the right Trach?* 2014. Available at: <http://bit.ly/ontherighttrach>. Accessed August 22, 2017.
2. Mallett J, Albarran JW, Richardson A, eds. *Critical Care Manual of Clinical Procedures and Competencies*. Oxford, UK: Wiley Blackwell; 2013.

Contents of unopened, undamaged package are clean, ready to use and for single-use. Refer to package insert for current warnings, indications, contraindications, precautions, and Instructions for Use.

Teleflex, the Teleflex logo, CrystalClear, PercuTwist, and Rüsch are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative.

Revised: 05/2018 © 2018 Teleflex Incorporated. All rights reserved. MC-004463 LA EN