



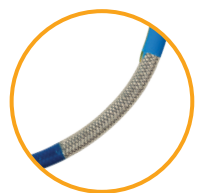

Pronto[®] V3 Extraction Catheter

Deliverability in Design

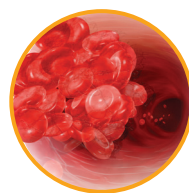

Deliverability and Vessel Protection

The Pronto® V3 Extraction Catheter was specially designed for deliverability during percutaneous thrombus aspiration with a uniquely engineered tip to help prevent vessel wall adhesion. The full-length braided shaft is engineered to enhance pushability and kink resistance while maintaining flexibility. This combination offers impressive deliverability when and where it's needed most.

Patented Silva tip

- Designed to prevent vessel wall adhesion
- Tapered distal tip designed to prevent vessel damage and improve navigation

PTFE lined wire lumen
for smooth guidewire movement

Shaft Construction

Proximal Shaft

Full-Length Braided Shaft


Distal Shaft

Hydrophilic coating
on the distal segment
promotes smooth vessel
navigation and trackability

Designed for Rapid Thrombus Removal

Uniformly Large Extraction Lumen

Engineered to maximize thrombus extraction force

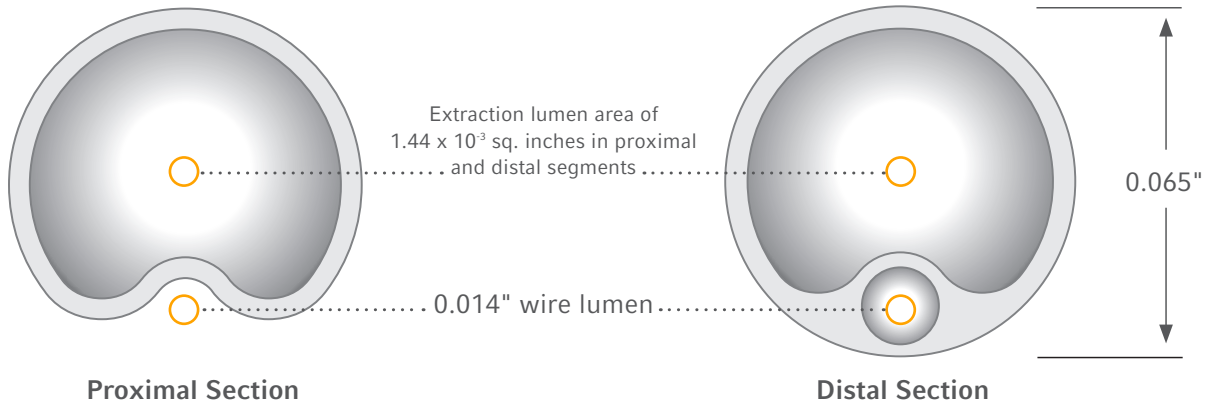


Aspiration Speed vs. Competition*


PRONTO® V3 CATHETER	EXPORT® ADVANCE	PRIORITYONE™ AC	EXPORT® AP
1.78 cc/s	1.73 cc/s	1.64 cc/s	1.55 cc/s

Tapered, Silva Tip vs. Open Tip

Designed to improve navigation and prevent vessel wall adhesion

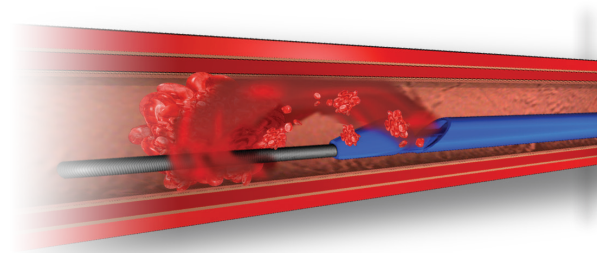

Patented Silva Tip⁺

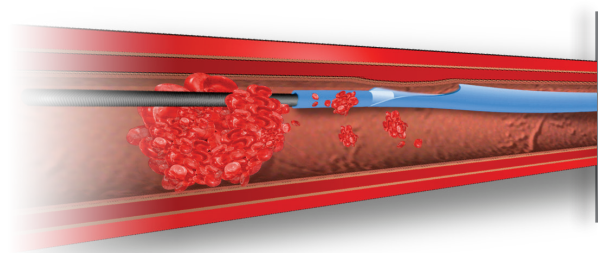

Standard Open-End Extraction Catheter Tip⁺

Specifications and Ordering Information

MODEL	DESCRIPTION	GUIDE CATHETER COMPATIBILITY	GUIDEWIRE COMPATIBILITY	WORKING LENGTH	RAPID EXCHANGE SECTION	BOX QUANTITY
5003	Pronto® V3 Extraction Catheter	≥6F (0.070" min. I.D.)	≤0.014"	140 cm	20 cm	1

Packaged in quantities of 1 unit per box.

The Pronto® V3 Extraction Catheter is indicated for the removal of fresh, soft emboli and thrombi from vessels in the coronary and peripheral vasculature.

*Testing completed by Teleflex. Data on file. Test results based on average of 10 aspiration cycles at an average extraction rate of 15 cc/s. Bench test results are not necessarily indicative of clinical performance.

+For illustrative purposes only; based on bench testing. Data on file. Bench test results are not necessarily indicative of clinical performance.

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow®, Deknatel®, Hudson RCI®, LMA®, Pilling®, Rüschi®, UroLift®, and Weck® – trusted brands united by a common sense of purpose.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Australia 1300 360 226

Austria +43 (0)1 402 47 72

Belgium +32 (0)2 333 24 60

Canada +1 (0)800 387 9699

China (Shanghai) +86 (0)21 6163 0965

China (Beijing) +86 (0)10 6418 5699

Czech Republic +420 (0)495 759 111

France +33 (0)5 62 18 79 40

Germany +49 (0)711 2090 8000

Greece +30 210 67 77 717

India +91 (0)44 2836 5040

Italy +39 0362 58911

Japan +81 (0)3 6632 3600

Korea +82 2 536 7550

Mexico +52 55 5002 3500

Netherlands +31 (0)88 00 215 00

New Zealand 0800 601 100

Poland +48 22 4624032

Portugal +351 22 541 90 85

Singapore (SEA non-direct sales countries) +65 6439 3000

Slovak Republic +421 (0)3377 254 28

South Africa +27 (0)11 807 4887

Spain +34 918 300 451

Switzerland +41 (0)31 818 40 90

United Kingdom +44 (0)1494 53 27 61

For more information, please visit teleflex.com.

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

Teleflex, the Teleflex logo, and Pronto are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. All other trademarks are trademarks of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 07/2019

© 2019 Teleflex Incorporated. All rights reserved. MC-005605