

Pilling™

The Pilling™ Vitality Line of Laryngoscopes



Teleflex™
Empowering the future of healthcare




Features and Benefits

- Modular accessories
- Color-coding
- All scopes laser compatible
- All scopes have jet ventilation port
- 3D printed solid stainless steel construction

Our most popular Dedo

The Dedo-Pilling Micro laryngoscopes are designed to provide the ability to examine the vocal cords at the anterior commissure in the majority of patients. The dimensions of the Dedo-Pilling laryngoscopes are narrow enough to permit moving the instrument off the incisors and down between the molars to permit visualization of the anterior commissure.

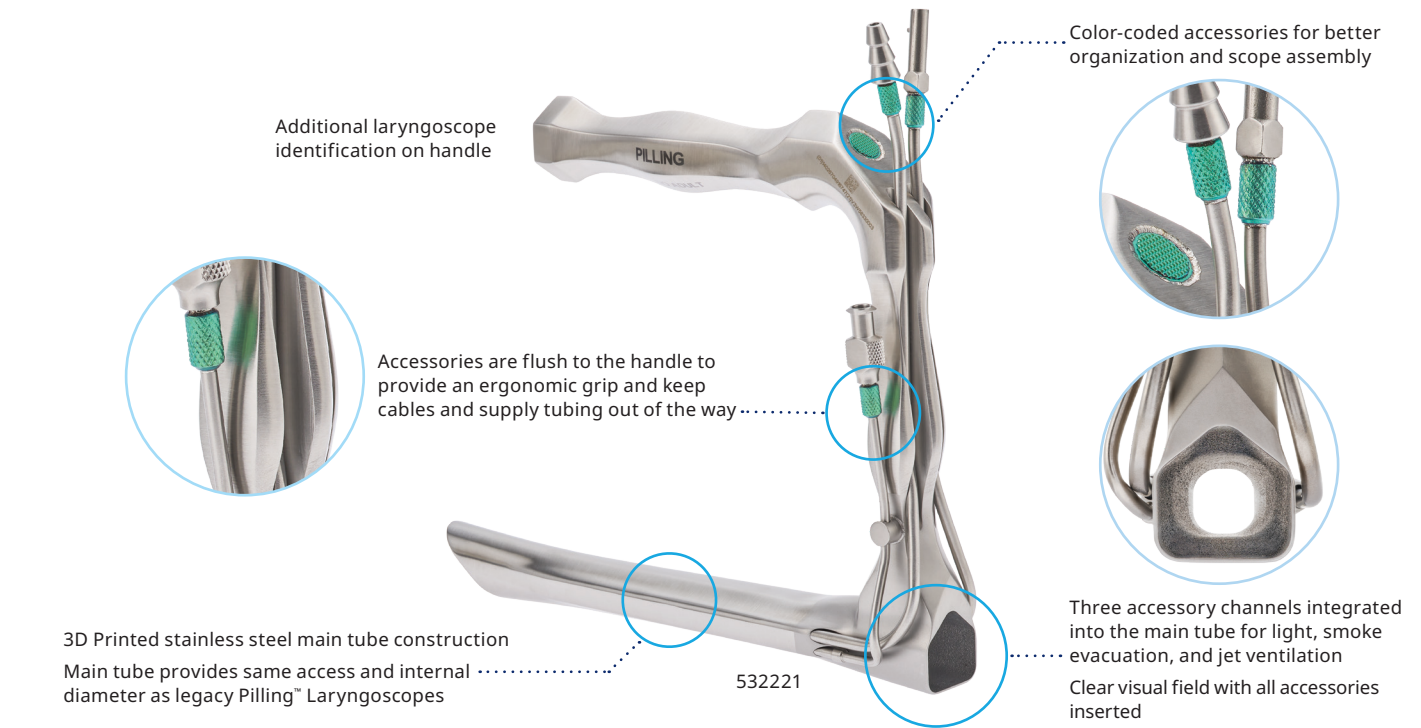

Pilling™ Vitality Dedo Micro Laryngoscopes

ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	INDIVIDUAL COMPONENTS				
					SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532221	Adult Dedo Laser Micro Laryngoscope	178 mm	10.7 mm	19.0 mm x 19.0 mm	532221S	537525	538525	539525	



532222



532224

532222	Child Dedo Micro Laser Laryngoscope	119 mm	7.1 mm	15.0 mm x 16.0 mm	532222S	537524	538524	539524	
532224	Adolescent Dedo Micro Laser Laryngoscope	145 mm	9.1 mm	18.0 mm x 17.0 mm	532224S	537524	538524	539524	

Ossoff-Karlan Family

The Vitality Ossoff-Karlan provide the same larger proximal width and height as legacy Pilling laryngoscopes, however the elimination of the two large oval light carrier channels reduce the outer profile by 15%.



Pilling™ Vitality Ossoff-Karlan Laryngoscopes


ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	INDIVIDUAL COMPONENTS				
					SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532231	Adult Ossoff-Karlan Laser Laryngoscope	178 mm	11.9 mm	22.0 mm x 20.0 mm	532231S	537525	538525	539525	
532237	Child Ossoff-Karlan Laser Laryngoscope	130 mm	9.8 mm	18.0 mm x 18.0 mm	532237S	537524	538524	539524	
532238	Adolescent Ossoff-Karlan Laser Laryngoscope	152 mm	9.7 mm	21.0 mm x 20.0 mm	532238S	537524	538524	539524	



532233	Ossoff-Karlan Posterior Commissure Laryngoscope	177 mm	12.6 mm	22.0 mm x 19.0 mm	532233S	537525	538525	539525	
532243	Dedo Laser Extra-Long Laryngoscope	197 mm	12.6 mm	20.0 mm x 20.0 mm	532243S	537525	538525	539525	


Ossoff-Pilling Micro-Laryngoscopes

Available in two proximal widths of 19 mm and 25 mm, all Ossoff-Pilling laryngoscopes provide the surgeon binocular vision while allowing for ease of introduction with a Holinger distal end. The Garrett-Ossoff-Pilling laryngoscope provides an additional 2 cm of reach and is designed for patients with a greater distance between their teeth and vocal cords.

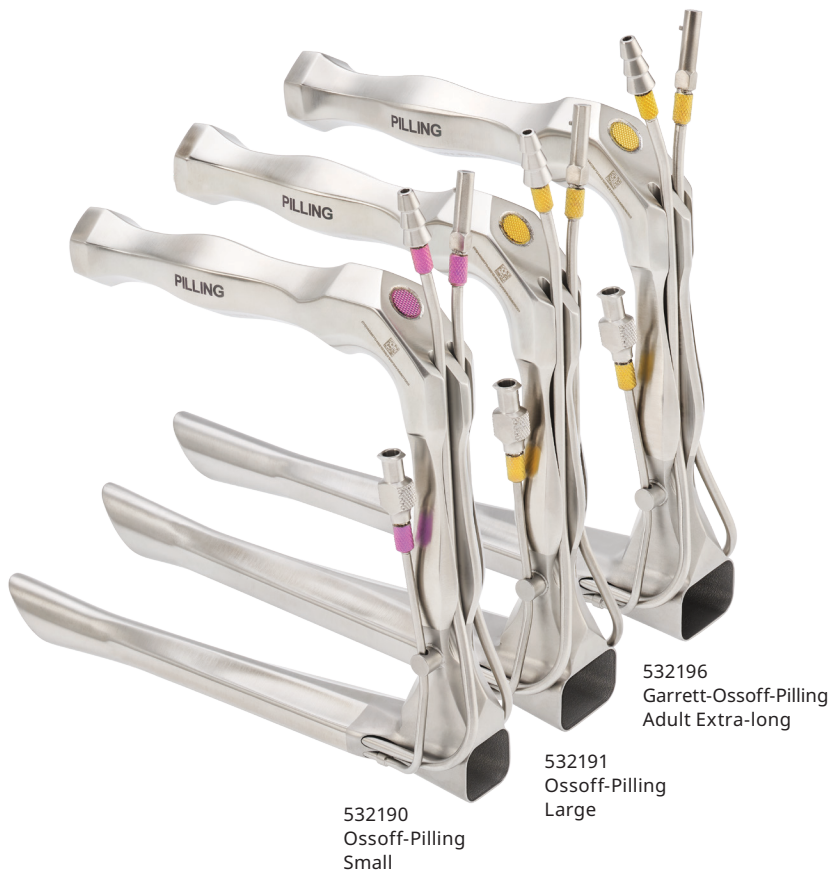

Upgrade your anterior scope to Ossoff-Pilling

Same anterior access, laser compatible scope that allows you to have binocular vision

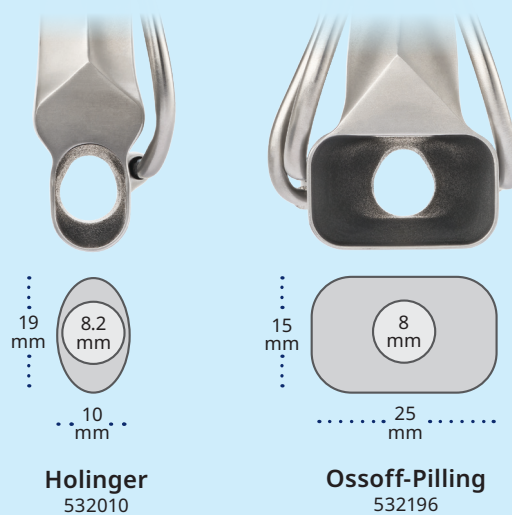


Pilling™ Vitality Ossoff-Pilling Laryngoscopes


ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	INDIVIDUAL COMPONENTS				
					SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532190	Small Ossoff-Pilling Laser Micro Laryngoscope	168 mm	7.0 mm	19.0 mm x 15.0 mm	532190S	537527	538527	539527	
532191	Large Ossoff-Pilling Laser Micro Laryngoscope	168 mm	8.0 mm	25.0 mm x 15.0 mm	532191S	537865	538865	539865	
532196	Extra-long Garrett-Ossoff-Pilling Laser Micro Laryngoscope	188 mm	9.0 mm	25.0 mm x 15.0 mm	532196S	537865	538865	539865	

Pilling™ Lindholm-Style Laser Laryngoscope

The Pilling™ Lindholm-Style Laryngoscope is the newest addition to the Pilling™ laryngoscope offering. All accessories are modular, allowing the surgeon to tailor the laryngoscope to the patient and condition being treated. With accessory integration in the main tube the Lindholm-Style Laryngoscope provides visualization of the posterior commissure, epiglottis, and cricoid region.

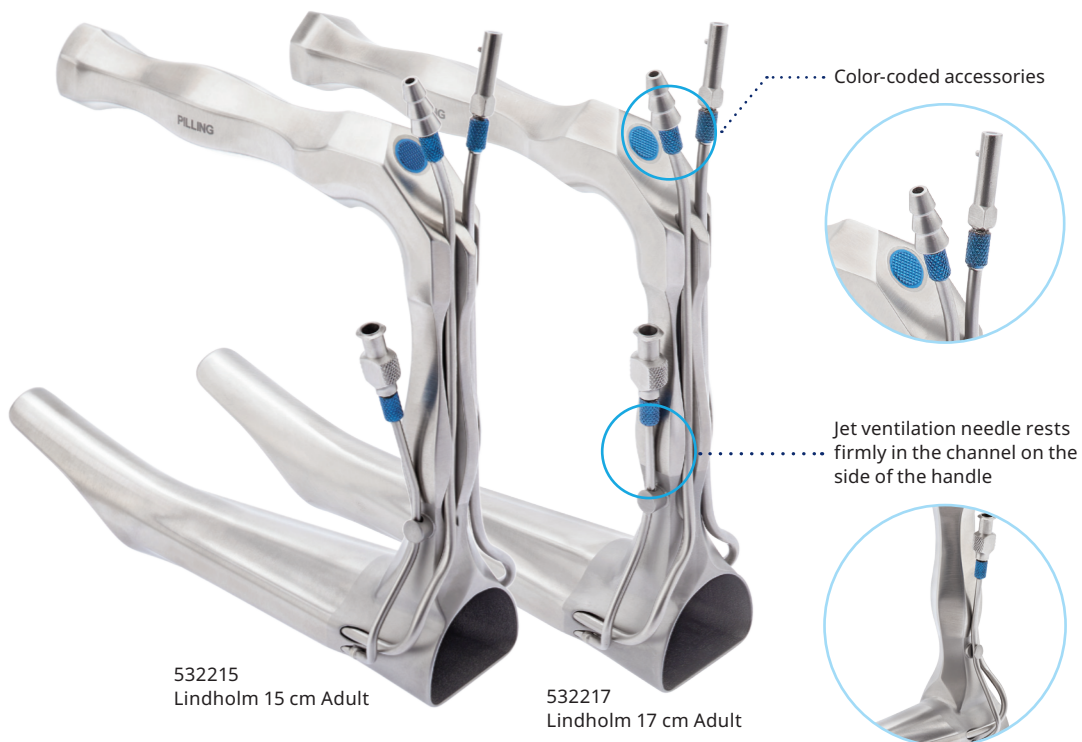

Minimal Distal Opening and Proximal Dimensions

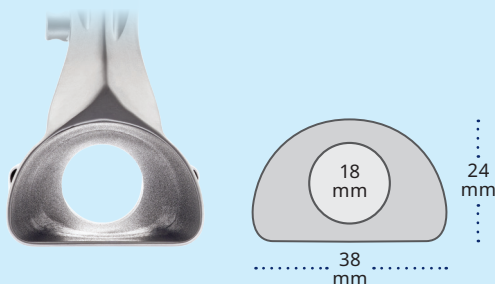

Pilling™ Vitality Lindholm-Style Laryngoscopes

ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	INDIVIDUAL COMPONENTS				
					SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532215	Adult Lindholm-Style Laser Laryngoscope 15 cm	150 mm	18 mm	38 mm x 24 mm	532215S	537215	538215	539215	
532217	Adult Lindholm-Style Laser Laryngoscope 17 cm	170 mm	18 mm	38 mm x 24 mm	532217S	537215	538215	539215	

Jako Laryngoscopes

Jako Micro-Laryngoscopes feature a wide proximal opening and a cutout at the distal end.



532206

Pilling™ Vitality Jako Laryngoscopes

ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	INDIVIDUAL COMPONENTS				
					SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532206	Adult Jako Laser Laryngoscope	178 mm	13.4 mm	29.0 mm x 20.0 mm	532206S	537865	538865	539865	Yellow
532209	Adolescent Jako-Cherry Laser Laryngoscope	149 mm	11.7 mm	25.0 mm x 13.0 mm	532209S	537527	538527	539527	Purple
532216	Child Jako Micro Laser Laryngoscope	120 mm	9.1 mm	20.0 mm x 13.0 mm	532216S	537524	538524	539524	Pink

Holinger Laryngoscope

The Holinger Anterior Commissure Laryngoscope incorporates the principle of the Yankauer post-nasal speculum. The distal end of the superior surface of the tube rises obliquely and flares laterally and the proximal portion drops obliquely. This increases the lumen and gives a larger field for work by elevating the tuberculum of the epiglottis.




532010

Pilling™ Vitality Holinger Laryngoscope

					INDIVIDUAL COMPONENTS				
ITEM NO.	DESCRIPTION	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532010	Holinger Adult Laryngoscope	170 mm	8.2 mm	10 mm x 14 mm	532010S	537028			

Pilling™ Vitality Next Gen Laryngoscopes Product Catalog

						INDIVIDUAL COMPONENTS				
ITEM NO.	DESCRIPTION	LEGACY ITEM NO.	OVERALL LENGTH	MIN I.D.	NOMINAL PROXIMAL DIMENSIONS	SCOPE ONLY	LIGHT CARRIER	SMOKE EVACUATION TUBE	JET VENTILATION	COLOR CODE
532010	Holinger Adult Laryngoscope	522010	170 mm	8.2 mm	10 mm x 14 mm	532010S	537028			
532190	Small Ossoff-Pilling Laser Micro Laryngoscope	522190	168 mm	7.0 mm	19.0 mm x 15.0 mm	532190S	537527	538527	539527	
532191	Large Ossoff-Pilling Laser Micro Laryngoscope	522191	168 mm	8.0 mm	25.0 mm x 15.0 mm	532191S	537865	538865	539865	
532196	Extra-long Garrett-Ossoff-Pilling Laser Micro Laryngoscope	522196	188 mm	9.0 mm	25.0 mm x 15.0 mm	532196S	537865	538865	539865	
532206	Adult Jako Laser Laryngoscope	522206	178 mm	13.4 mm	29.0 mm x 20.0 mm	532206S	537865	538865	539865	
532209	Adolescent Jako-Cherry Laser Laryngoscope	522209	149 mm	11.7 mm	25.0 mm x 13.0 mm	532209S	537527	538527	539527	
532215	Adult Lindholm-Style Laser Laryngoscope 15 cm		150 mm	18 mm	38 mm x 24 mm	532215S	537215	538215	539215	
532216	Child Jako Micro Laser Laryngoscope	522216	120 mm	9.1 mm	20.0 mm x 13.0 mm	532216S	537524	538524	539524	
532217	Adult Lindholm-Style Laser Laryngoscope 17 cm		170 mm	18 mm	38 mm x 24 mm	532217S	537215	538215	539215	
532221	Adult Dedo Laser Micro Laryngoscope	522221, 522225P, 522228	178 mm	10.7 mm	19.0 mm x 19.0 mm	532221S	537525	538525	539525	
532222	Child Dedo Micro Laser Laryngoscope	522222	119 mm	7.1 mm	15.0 mm x 16.0 mm	532222S	537524	538524	539524	
532224	Adolescent Dedo Micro Laser Laryngoscope	522224	145 mm	9.1 mm	18.0 mm x 17.0 mm	532224S	537524	538524	539524	
532231	Adult Ossoff-Karlan Laser Laryngoscope	522231	178 mm	11.9 mm	22.0 mm x 20.0 mm	532231S	537525	538525	539525	
532233	Ossoff-Karlan Posterior Commissure Laryngoscope	522233	177 mm	12.6 mm	22.0 mm x 19.0 mm	532233S	537525	538525	539525	
532237	Child Ossoff-Karlan Laser Laryngoscope	522237	130 mm	9.8 mm	18.0 mm x 18.0 mm	532237S	537524	538524	539524	
532238	Adolescent Ossoff-Karlan Laser Laryngoscope	522238	152 mm	9.7 mm	21.0 mm x 20.0 mm	532238S	537524	538524	539524	
532243	Dedo Laser Extra-Long Laryngoscope	522243	197 mm	12.6 mm	20.0 mm x 20.0 mm	532243S	537525	538525	539525	

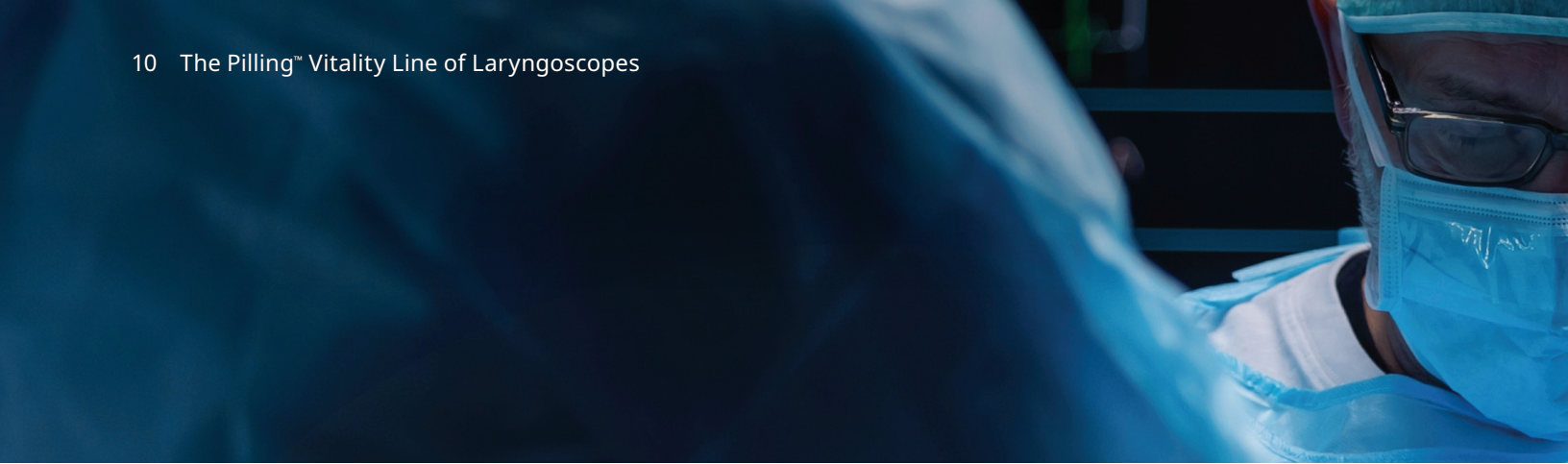

Additional Accessories

Laryngoscope Holders / Table Mounted Trays

ITEM NO.	DESCRIPTION
502245	Lewy-style Laryngoscope Holder
502248	Tobey Laryngoscope Holder with Quick Release
502260	Boston University Suspension System
502261	Boston University Laryngoscope Handle Extension
502265	Dedo Table Mounted Mayo Tray



Tobey Laryngoscope Holder with Quick Release
502248

Lightsource and Accessories

Headlight

ITEM NO.	DESCRIPTION
528339PS	Variable 20-110 mm Spot with Bifurcated 9' Fused Fiber Cable

Lightsources

ITEM NO.	DESCRIPTION
528300PS	Xenon 300 Watt Digital Lightsource, 5 Port Turret

Floorstand

ITEM NO.	DESCRIPTION
528304PS	Floorstand, Including Base/Pole/Wheel Base



300 Watt Digital Xenon Lightsource



Fiberoptic Cable Compatibility Guide

LIGHTSOURCE ENDFITTING	INSTRUMENT ENDFITTINGS				
	ACMI® FEMALE INSTRUMENT	BIFURCATED (2) PILLING™ ROTATING	PILLING™ ROTATING BAYONET (PILLING BLUE)		STORZ® THREADED INSTRUMENT
DESCRIPTION	8' LENGTH	8' LENGTH	6' LENGTH	8' LENGTH	8' LENGTH
Pilling™ Lightsource Endfitting (P)			P356PRA	P358PRA	
ACMI® Lightsource Endfitting (A)	A358AFA			A358PRA	A358STA A508STA*
Olympus® Lightsource Endfitting				O358PRA	
Storz® Lightsource Endfitting (S)				S358PRA	
Wolf® Lightsource Endfitting (W)		W358PXA		W358PRA	

Standard Bundle Size is 3.5 mm.
* 5.0 mm Bundle Size

Fiberoptic Cable

- Step 1: Select lightsource endfitting
- Step 2: Select bundle size
- Step 3: Select length
- Step 4: Select instrument endfitting



To view Pilling™ Vitality Laryngoscope resources go to: myTeleflex.com

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

Teleflex, the Teleflex logo, and Pilling are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

All other trademarks are trademarks of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 07/2024.

© 2024 Teleflex Incorporated. All rights reserved.

MC-009967 · REV 0