

The Chocolate XD[®] PTCA Balloon Catheter

Now Available from Teleflex

Designed for Predictable, Uniform, Atraumatic Dilatation

The Chocolate XD[®] PTCA Balloon Catheter is a specialty angioplasty balloon, used in the pre-dilation and treatment of coronary lesions. The balloon's proprietary nitinol constraining structure creates "pillows" and "grooves" that are designed to provide predictable, uniform, and atraumatic dilatation. The "pillows" are designed to provide predictable vessel dilatation without cutting or scoring, and the "grooves" allow for stress relief, potentially minimizing the traumatic angioplasty effect.



Features

Designed for vessel preparation

Expanded angioplasty treatment options

Designed to deliver controlled dilatation while minimizing vessel wall trauma

Multiple sizes to accommodate a wide range of complex lesions

The Chocolate XD[®] PTCA Balloon Catheter

Now Available from Teleflex


Indications for Use

The Chocolate XD[®] PTCA Balloon Dilatation Catheter is indicated for balloon dilatation of the stenotic portion of coronary artery or bypass graft stenosis for the purpose of improving myocardial perfusion including in-stent restenosis.

Ordering Information

CATALOG NUMBER	GUIDEWIRE COMPATIBILITY	CATHETER LENGTH	BALLOON DIAMETER	BALLOON LENGTH	NOMINAL PRESSURE	RATED BURST PRESSURE	FRENCH SIZE COMPATIBILITY
CX14-140-2010 RX	0.014"	140 cm	2.0 mm	10 mm	10 atm	16 atm	5 Fr
CX14-140-2015 RX	0.014"	140 cm	2.0 mm	15 mm	10 atm	16 atm	5 Fr
CX14-140-2020 RX	0.014"	140 cm	2.0 mm	20 mm	10 atm	16 atm	5 Fr
CX14-140-2510 RX	0.014"	140 cm	2.5 mm	10 mm	10 atm	16 atm	5 Fr
CX14-140-2515 RX	0.014"	140 cm	2.5 mm	15 mm	10 atm	16 atm	5 Fr
CX14-140-2520 RX	0.014"	140 cm	2.5 mm	20 mm	10 atm	16 atm	5 Fr
CX14-140-3010 RX	0.014"	140 cm	3.0 mm	10 mm	10 atm	16 atm	6 Fr
CX14-140-3015 RX	0.014"	140 cm	3.0 mm	15 mm	10 atm	16 atm	6 Fr
CX14-140-3020 RX	0.014"	140 cm	3.0 mm	20 mm	10 atm	16 atm	6 Fr
CX14-140-3510 RX	0.014"	140 cm	3.5 mm	10 mm	10 atm	16 atm	6 Fr
CX14-140-3515 RX	0.014"	140 cm	3.5 mm	15 mm	10 atm	16 atm	6 Fr
CX14-140-3520 RX	0.014"	140 cm	3.5 mm	20 mm	10 atm	16 atm	6 Fr

Packaged in quantities of one unit per box.

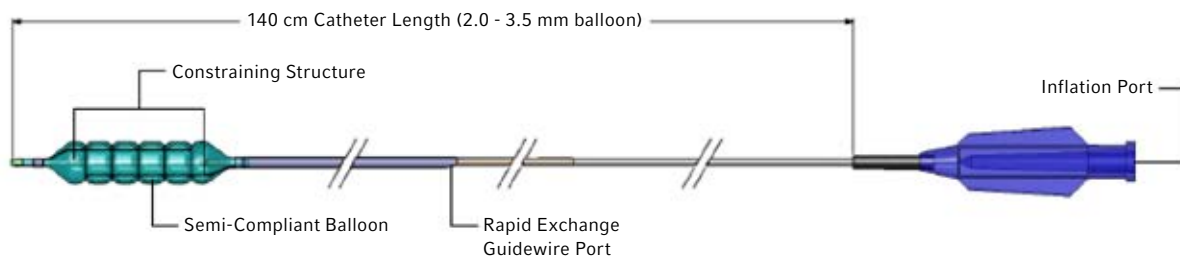

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions. CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

The Chocolate XD[®] PTCA Balloon Dilatation Catheter is contraindicated for unprotected left main coronary artery, coronary artery spasm in the absence of significant stenosis and crossing through stent struts.

Chocolate XD is a registered trademark of TriReme Medical LLC.

Teleflex and the Teleflex logo are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 06/2018.

© 2018 Teleflex Incorporated. All rights reserved. MC-004616 Rev 0